



## CLARUS 700 from ZEISS


HD Ultra-widefield Fundus Imaging with Fluorescein Angiography



Seeing beyond

# Expanding insights with ultra-wide imaging.

## ZEISS CLARUS 700

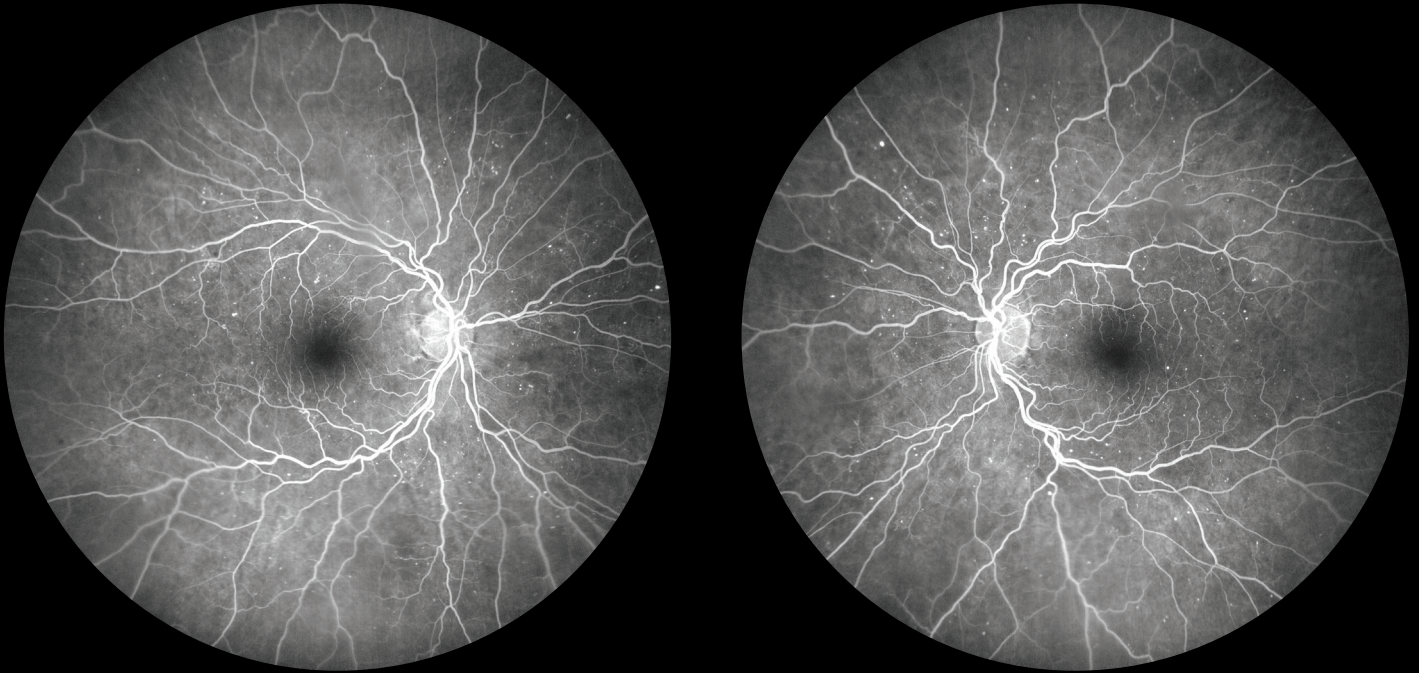


Fluorescein angiograms with non-proliferative diabetic retinopathy, illustrating localized dilations of retinal capillaries (microaneurysms) and areas of peripheral non-perfusion.

## Unsurpassed image quality with fluorescein angiography.

CLARUS® 700 from ZEISS allows you to capture clear and accurate images from the macula to the far periphery, all with a single instrument that combines:

- Ultra-wide field of view
- True Color imaging from broad spectrum LED scans
- Exceptional resolution
- Fluorescein Angiography (FA)
- Advanced imaging features

ZEISS CLARUS 700 is a truly comprehensive imaging system developed for eye care specialists, helping deliver state-of-the-art care to their patients.

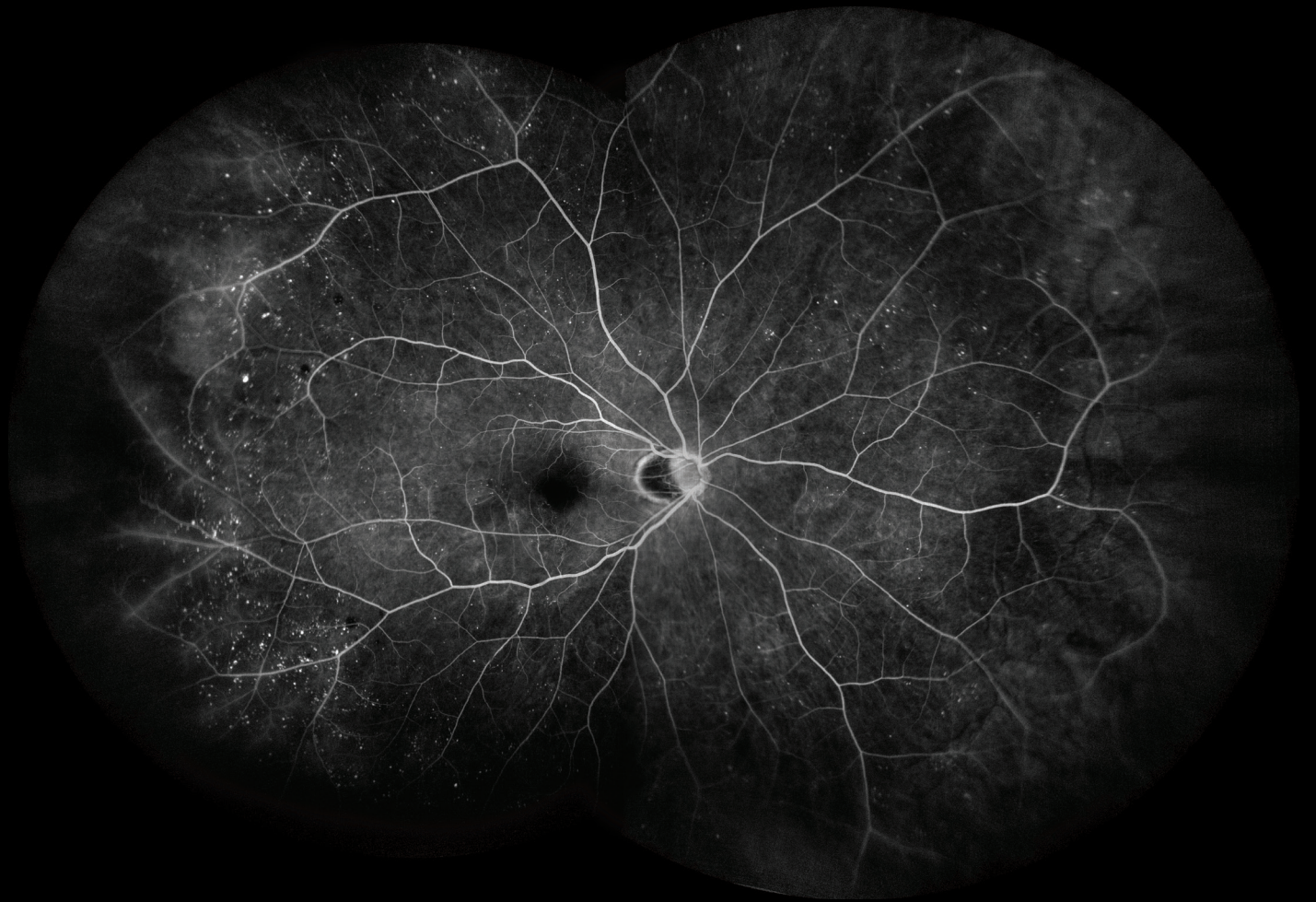

Montage fluorescein angiography with non-proliferative diabetic retinopathy, presenting the finest details at the foveal avascular zone and offering an exceptional rendering of the smallest microaneurysms across the image—from the fovea to the periphery.

**COLOR**

Capture True Color to assist with differential diagnosis.

**CLARITY**

See high-resolution details from the posterior pole to the periphery.

**COMPLETE**


Comprehensive in every way to maximize workflow efficiency.



# Now, ultra-wide fundus imaging with True Color and unmatched clarity. In one complete system.

## True Color Imaging

Powered by **Broad Line Technology**, the ZEISS CLARUS 700 captures images that closely resemble the coloration of the fundus as seen during clinical examination.

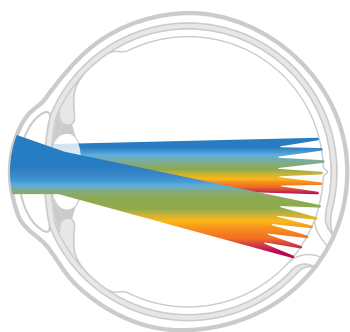

*Unlike CSLO (confocal scanning laser), Broad Line Technology enables the combination of ultra-wide fields of view and a full range of retinal imaging modes to generate images with high dynamic range, contrast, resolution and natural colors through sequential illumination of broad-spectrum red, green and blue light emitting diodes.<sup>1</sup>*

<sup>1</sup> Data on file.

## A Comprehensive Imaging System

Now you can manage all fundus imaging modalities without compromising on clarity—viewing high resolution in ultra-widefield.

- Image from the superior and inferior retina with less peripheral distortion
- Capture clear detail of vessel structure from early to late phase of fluorescein angiography
- **AutoBright** control automatically optimizes the angiogram series preserving change in signal


Combining ultra-widefield imaging with True Color, excellent clarity and a full suite of imaging modalities, ZEISS CLARUS 700 empowers you with features and capabilities that maximize workflow efficiency.

- Quickly and easily compare images over time and between image capture modes
- Provide a comfortable patient experience that ensures image integrity, with ergonomic chin and head rests to swivel motion and live IR preview


## Complete suite of imaging modalities




True Color with RGB Channel Separation

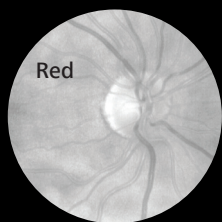


Red channel: reveals the choroid in more detail. This may be helpful in visualizing choroidal lesions such as nevi or tumors.

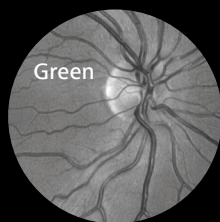

Green channel: provides excellent contrast of the retina, especially of vasculature and hemorrhages.

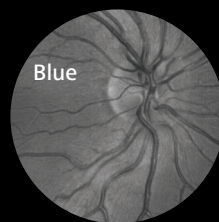


Blue channel: increases visibility of anterior retinal layers, allowing easier visualization of retinal nerve fiber layers defects and epiretinal membranes.

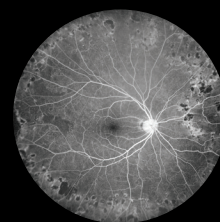


Fluorescein Angiography of proliferative diabetic retinopathy




FAF-Green image of dry age-related macular degeneration

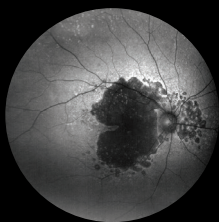

FAF-Blue image of geographic atrophy

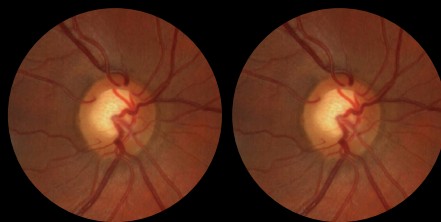

Stereo image pairs can be captured for stereoscopic evaluation of the fundus.

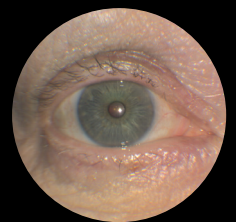



External Eye

## True advancement in disease management.

By allowing you to visualize to the far periphery—and in multiple imaging modalities—ZEISS CLARUS 700 can document indications of ocular disease that occur in various regions of the eye and present differently depending on the imaging modality.

### Proliferative Diabetic Retinopathy


Early phase fluorescein angiogram: Visualize macular ischemia, capillary nonperfusion and intraretinal microvascular abnormalities in excellent detail with high-resolution imaging.

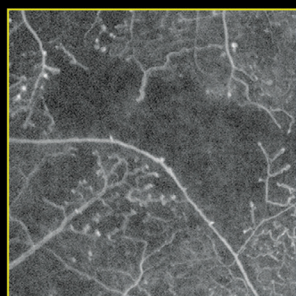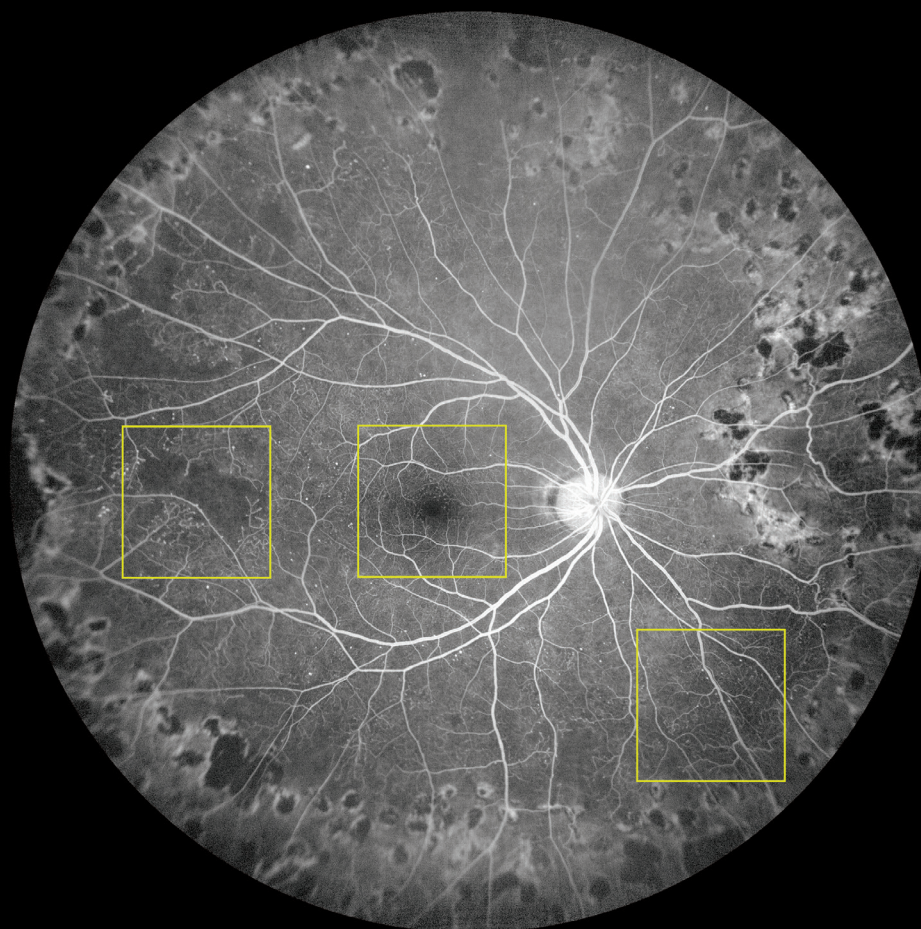

*Capillary nonperfusion*




*Macular ischemia*

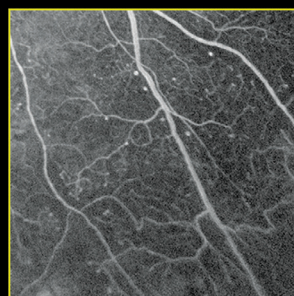

*Intraretinal microvascular abnormalities*

## Macular Telangiectasia



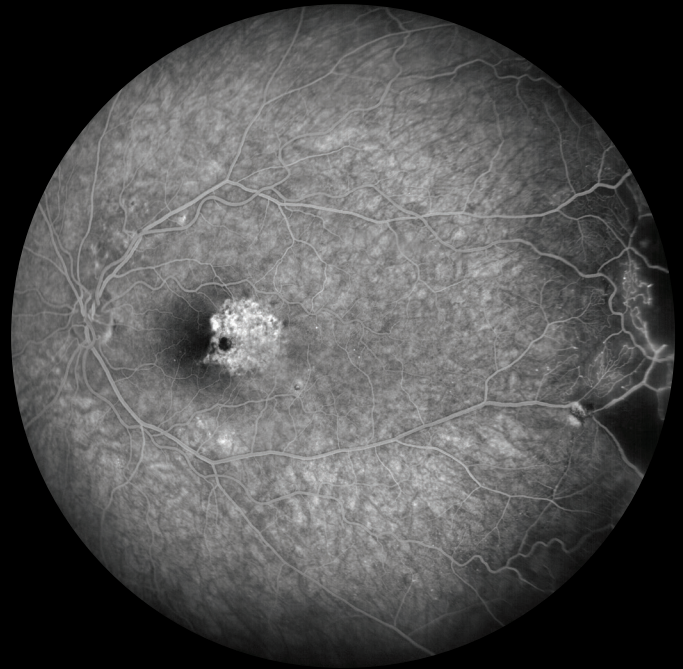


Mid-phase fluorescein angiography image of an eye with macular telangiectasia. Wide-field fluorescein angiography captures leakage in the macula, its associated microaneurysms and non-perfusion in the far temporal periphery.

## Dry AMD



FAF-Green image of an eye with central geographic atrophy in advanced Dry AMD, highlighting the loss of retinal pigment epithelium at the macula.

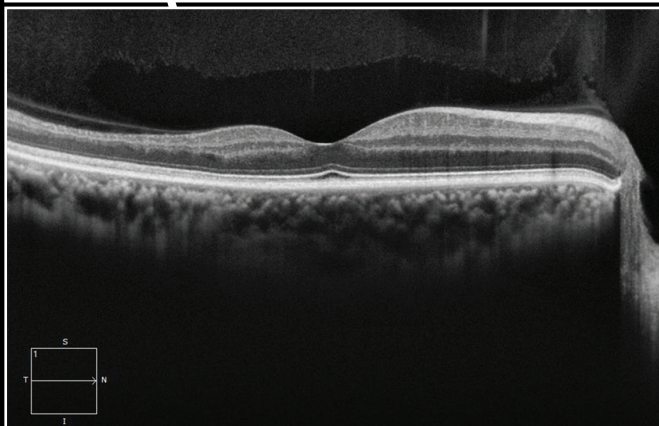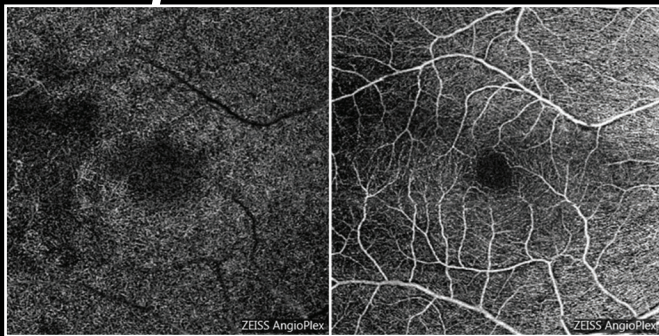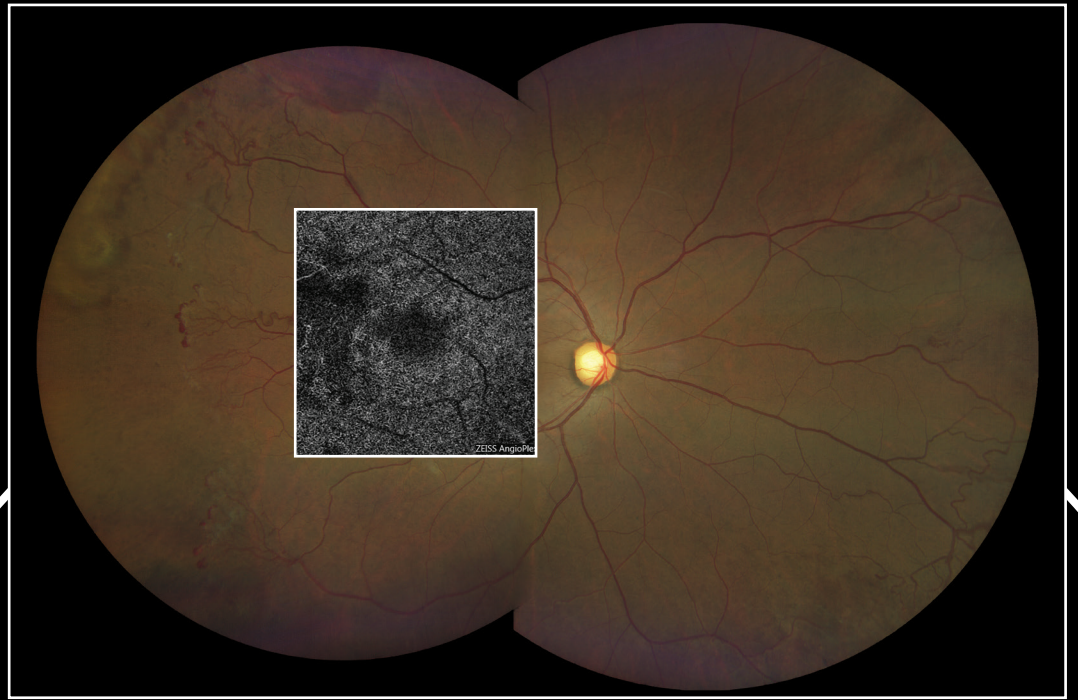

ZEISS CLARUS 700 ultra-wide fluorescein angiography shows you the extensive sea-fan neovascularization and retinal ischemia in the peripheral retina in an eye with proliferative sickle cell retinopathy.




# Integrated Diagnostic Imaging platform from ZEISS.


See the whole picture.

Key to meeting current challenges in eye care is the ability to capture, integrate and transform high-quality data into meaningful analyses that enhance practice workflow and improve patient care.

The ZEISS Integrated Diagnostic Imaging combines exam data from gold-standard devices like CLARUS ultra-widefield fundus imaging and CIRRUS™ HD-OCT from ZEISS and presents critical information from multiple sources into a single integrated point-of-view for more efficient and insightful treatment decisions.

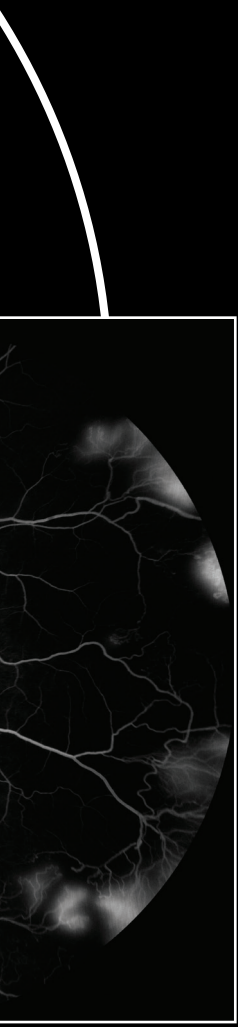

Retina Workplace

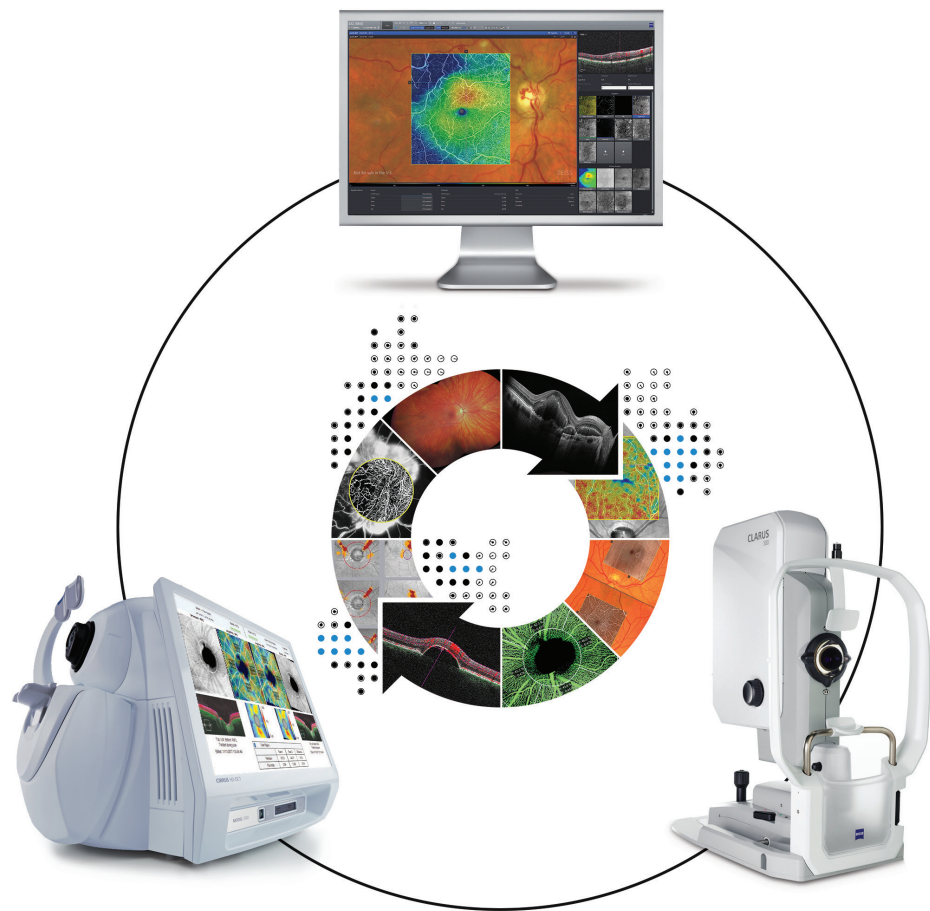


CIRRUS™ HD-OCT

CLARUS® 700

# Advanced features to help you capture your best images.


## PrecisionFocus

Quickly see the details in regions of interest by selecting where to optimize focus, without losing the macula focal point.

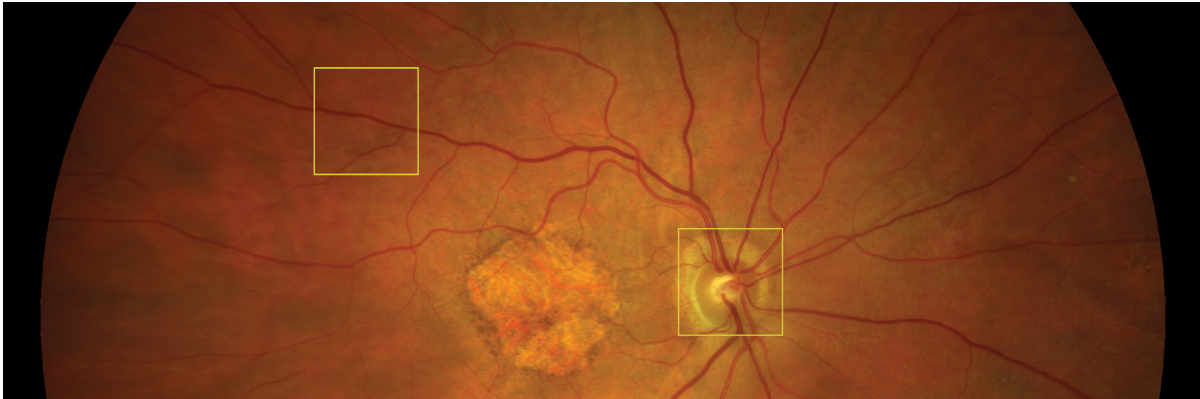


## AutoBright

Spend time analyzing images rather than adjusting them. ZEISS CLARUS 700 automatically optimizes the brightness of the image sequence throughout the angiogram, while still preserving the change in signal. And with the extremely large dynamic range, you'll never be at risk of saturating the image.

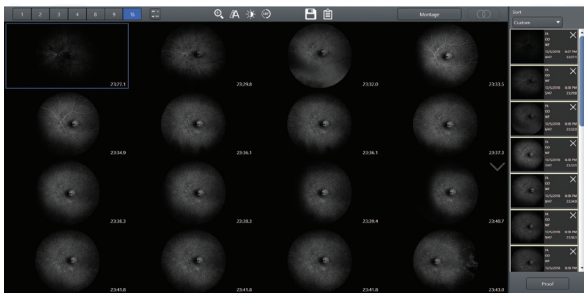


Original

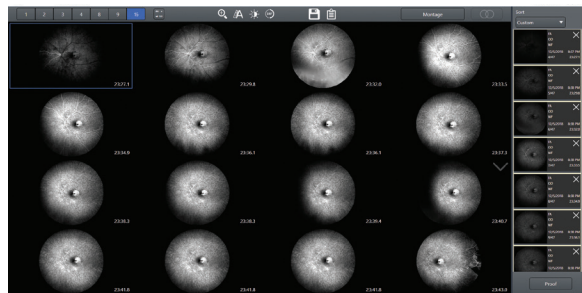

AutoBright

## GazePoint

Find the patient's gaze angle quickly and accurately. CLARUS 700 uses AI to automatically find the optic nerve head and accurately derive the patient's gaze rather than relying on internal fixation.

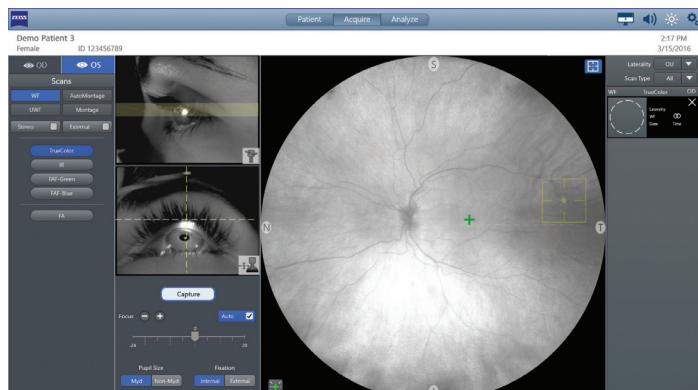

# Technical Specifications

## CLARUS 700 from ZEISS

### Parameters

#### Imaging Modes:

- True Color (with Red, Green and Blue channel separation)
- Fluorescein Angiography
- Autofluorescence-Green
- Autofluorescence-Blue
- External eye image (ocular surface)
- Stereo

#### Field of View (measured from the center of the eye):

- Widefield (one image) 133°
- Ultra-widefield (two images) 200°
- Montage (up to six images) up to 267°

#### Resolution:

- Optical 7.3 µm

**Minimum Pupil Diameter:** 2.5 mm

**Working Distance:** 25 mm (patient's eye to front lens)

**Compensation for ametropia:** - 24 D to + 20 D continuous

#### Light Sources:

- Red LED 585 - 640 nm
- Green LED 500 - 585 nm
- Blue LED 435 - 500 nm
- Infrared laser diode 785 nm

#### Automatic Operations:

- Auto-focus Auto Montage
- Auto-gain Auto-laterality

#### Aquisition Speed:

- Live IR Preview 10 frames/second
- Image Capture ≤ 0.2 seconds

### Instrument Specifications

**Acquisition Device Weight:** ~23.6 kg

**Acquisition Device Dimensions (W x D x H):** 362 mm x 546 mm x 676 mm

#### Instrument Table:

- Description Wheelchair accessible, electronic lift
- Table Dimensions (W x D x H) 916 mm x 615 mm x 711 - 925 mm
- Weight ~38 kg

#### Instrument Input Power:

- External Power Supply 100-240VAC, 50/60 Hz
- Electrical Class IEC 60601-1 Class I

### At-Instrument Computer

<b>Monitor:</b>	22" Full HD MVA LCD with LED Backlight	<b>Touch Screen:</b>	Capacitive, Multi-Touch
<b>Resolution:</b>	1920 x 1080	<b>RAM:</b>	32GB
<b>Processor:</b>	Intel® 6th Generation Core i5-6500TE	<b>Input/Output:</b>	USB 3.0 x 4; RS-232 x 2; 1.5 kV Isolated Gigabit Ethernet Port x 2; HDMI; and DisplayPort
<b>Hard Drive:</b>	2 TB (minimum 200,000 images)	<b>Operating System:</b>	Windows 10
<b>Dimensions (W x D x H):</b>	21.5" (54.6 cm) x 2.5" (6.4 cm) x 13.75" (34.9 cm)		
<b>Weight:</b>	~ 8.5 kg	<b>Mounting:</b>	VESA 75/100 mm



Download the ZEISS Image Library App directly from the App Store.  
Explore a wide selection of modalities such as ultra-widefield and OCTA.



0297  
CLARUS 700



**Carl Zeiss Meditec, Inc.**  
5300 Central Parkway  
Dublin, CA 94568  
USA  
[www.zeiss.com/us/clarus700](http://www.zeiss.com/us/clarus700)  
[www.zeiss.com/us/med](http://www.zeiss.com/us/med)



**Carl Zeiss Meditec AG**  
Goeschwitzer Strasse. 51-52  
07745 Jena  
GERMANY  
[www.zeiss.com/clarus700](http://www.zeiss.com/clarus700)  
[www.zeiss.com/med/contacts](http://www.zeiss.com/med/contacts)

en-INT\_31\_010\_0036111 CZ VI/2021

The contents of the brochure may differ from the current status of approval of the product or service offering in your country. Please contact our regional representatives for more information. Subject to changes in design and scope of delivery and due to ongoing technical development. CLARUS 700 and CIRRUS-HDD-OCT are either trademarks or registered trademarks of Carl Zeiss Meditec AG or other companies of the ZEISS Group in Germany and / or other countries.  
© Carl Zeiss Meditec, Inc., 2021. All rights reserved.